


# QUICK START GUIDE




## CHECK LIST

Ensure you have the following materials ready to run a test:

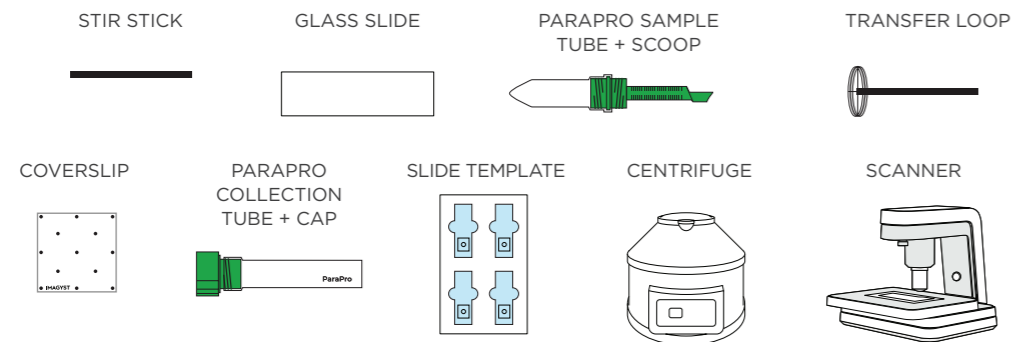


## 1 COLLECT FAECAL MATTER

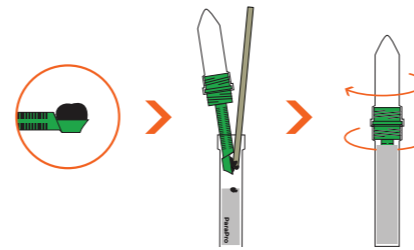


- Unscrew the **tube cap** from the **collection tube**
- Add 1 scoop, about 1 gram, of faecal matter to the **collection tube**
- Use a **stir stick** to remove the matter from the **scoop** and mix it
- Screw the **sample tube** onto the **collection tube**
- Continue twisting until the tubes are firmly screwed together

## 2 PREPARE FAECAL SAMPLE

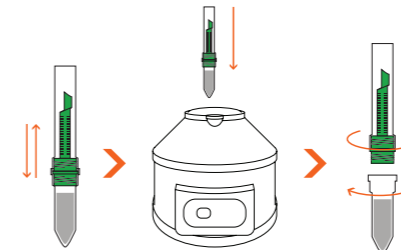


- Shake the tubes and then place them in a centrifuge with the collection tube **facing down** and the sample tube with scoop inside it **facing up**
- Spin the sample for **2 minutes** at 500 RCF
- Unscrew the **collection tube** with the **scoop** and discard both
- Place the **sample tube** in a tray

## 3 ADD SAMPLE TO SLIDE

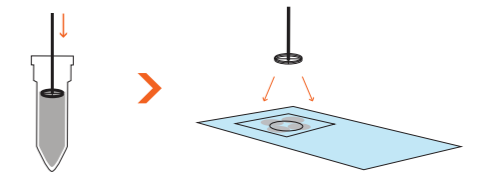


- Place a **glass slide** in the **slide template** with the name or barcode towards the top
- Pull material from the top of the **sample tube** using the **transfer loop**
- Place 2-4 loops of faecal matter onto the **glass slide** around the edge of the target circle

## 4 PLACE COVERSGLIP

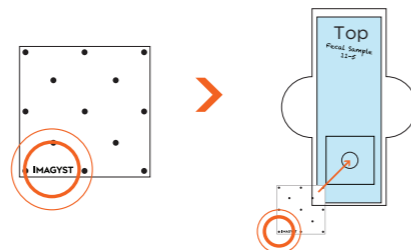

- Orient the **coverslip** with the **IMAGYST** text readable in the lower left-hand corner of the **coverslip**
- Place the coverslip on the **glass slide** over the targeted area of the **slide template**

## 5 INITIATE TEST

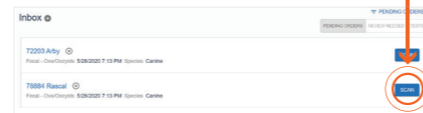


- Select **Scan** for the correct sample; Sample information will come pre-populated from any FUSE-connected software system
- If no integration is available, select the **Add Sample** icon (+), enter in the sample information, select **Create**, and then select **Scan**
- Select a **scanner**, and make sure the tray is open (select **Open Tray** if necessary)

## 6 LOAD AND SCAN SLIDE



- Unlock the slide holder and place the **glass slide** (with **coverslip**) on the slide tray. Place any label or writing towards the **right-side** of the slide tray and lock the slide holder
- Select **Scan Slide**

## 7 REVIEW SCAN

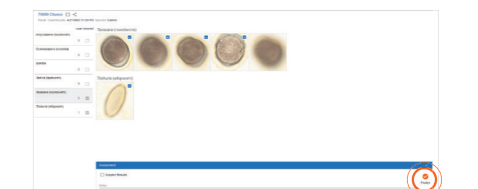

- Open the sample results, review any captured images, and select examples for inclusion in the report
- Add any comments and select **Finalise**
- Select the **View Report** icon to view a shareable PDF or review results in FUSE-connected software system

CONTACT ZOETIS TECHNICAL SUPPORT FOR ADDITIONAL INFORMATION

1800 270 727

DxSupport.au@zoetis.com

OR WATCH INSTRUCTIONAL VIDEOS AT:

[vetscanimagyst.com.au/resources](https://vetscanimagyst.com.au/resources)